



CANKADO
Achieving more together



Improving COVID-19 patient
care using the

CANKADO COVID-19 Caregiver Cockpit

Physician Instruction Manual

Version 1.0 | Date: 24.03.2020

CANKADO COVID-19 Caregiver Cockpit is registered as a Class I medical device within the European Union (registration number DE/CA59/BS 1599/2020-R/Hi) and is compliant with the FDA classification for Mobile Medical Devices (2015) Appendix B.

HCP Pro App Link



Content

Using the web application

01 Cockpit overview	3
----------------------------	----------

02 How to connect to a patient	7
---------------------------------------	----------

03 How to register	9
---------------------------	----------

Using the HCP Pro APP on a smartphone


01 App overview	14
------------------------	-----------

02 How to connect to a patient	19
---------------------------------------	-----------

03 How to register › go to web-registration	22
--	-----------

Using the web application

Cockpit overview

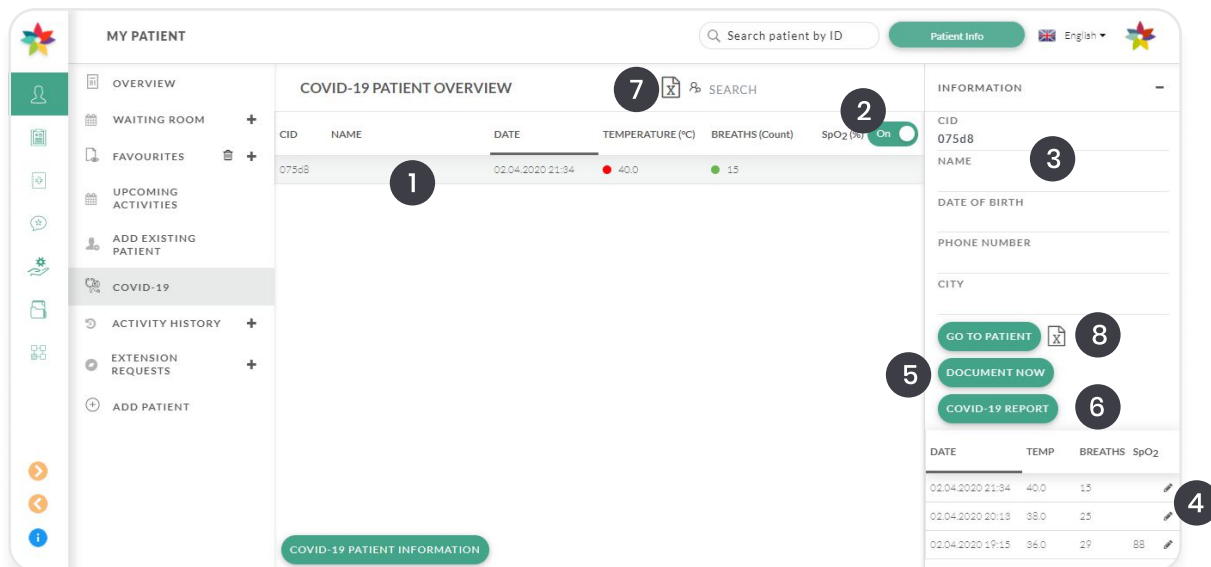

All COVID-19 related features are displayed in a separate cockpit. You have to select COVID-19 in the menu bar to reach here.

1: Patient list

All your patient will be shown here along with the parameters. by selecting each patient you can see related data in the right side column.

Colour code	Temperature	Breaths	SpO2
green	< 37.5 °C	< 19	> 94 %
Orange	37.5 – 38.5 °C	19 – 26	90 – 94 %
red	> 38.5 °C	> 26	< 90 %

What is the meaning of the arrows behind the values?

Arrows appear behind the values if there is a change compared to the value of yesterday. An oblique arrow indicates small changes and horizontal arrows significant changes.

If values stay unchanged or if no value for the day before is entered, there will be no arrow.

2: SpO₂ values on/off option for patients.

If you are monitoring patients who are not able to determine their SpO₂ value (e.g. because they are at home without a pulse oximeter), you can turn off SpO₂ documentation for all your patients.

3: Basic information about the selected patient

Here you can see basic patient information like CID (Cankado Patient ID), name, date of birth, phone number and city.

4: Quick edit of patient data

MY PATIENT

Search patient by ID

Patient Info English

OVERVIEW

WAITING ROOM +

FAVOURITES +

UPCOMING ACTIVITIES

ADD EXISTING PATIENT

COVID-19

075D8

ACTIVITY HISTORY +

EXTENSION REQUESTS +

ADD PATIENT

COVID 19

02.04.2020 21:34

SAVE

Item Name Value Unit Reference Range

Temperature 40 °C 35 - 42 °C

Breaths 15 count 12 - 30 count

SpO2 % 0 - 100 %

ITEM INFO

PID

CID 075d8

NAME

DATE OF BIRTH

PHONE NUMBER

CITY

DELETE

After selecting [4] in Cockpit Overview as depicted on page 3 you are able to edit the corresponding data entry.

- 1) Select the field you wish to edit and enter the new value.
- 2) Click on "Save" to continue.

5: Document COVID-19 status for a patient

MY PATIENT

Search patient by ID

Patient Info English

OVERVIEW

WAITING ROOM +

FAVOURITES +

UPCOMING ACTIVITIES

ADD EXISTING PATIENT

COVID-19

075D8

ACTIVITY HISTORY +

EXTENSION REQUESTS +

ADD PATIENT

COVID 19

04.04.2020 Select Time

SAVE

Item Name Value Unit Reference Range

Temperature °C 35 - 42 °C

Breaths count 12 - 30 count

SpO2 % 0 - 100 %

ITEM INFO

PID

CID 075d8

NAME

DATE OF BIRTH


PHONE NUMBER

CITY

After selecting [5] in Cockpit Overview as depicted on page 3 you can add new data for the corresponding patient.

- 1) Click into the appropriate fields and enter values. By default, the current date and time are used as the time of measurement.
- 3) Click on "Save" to continue.


6: Download COVID-19 report

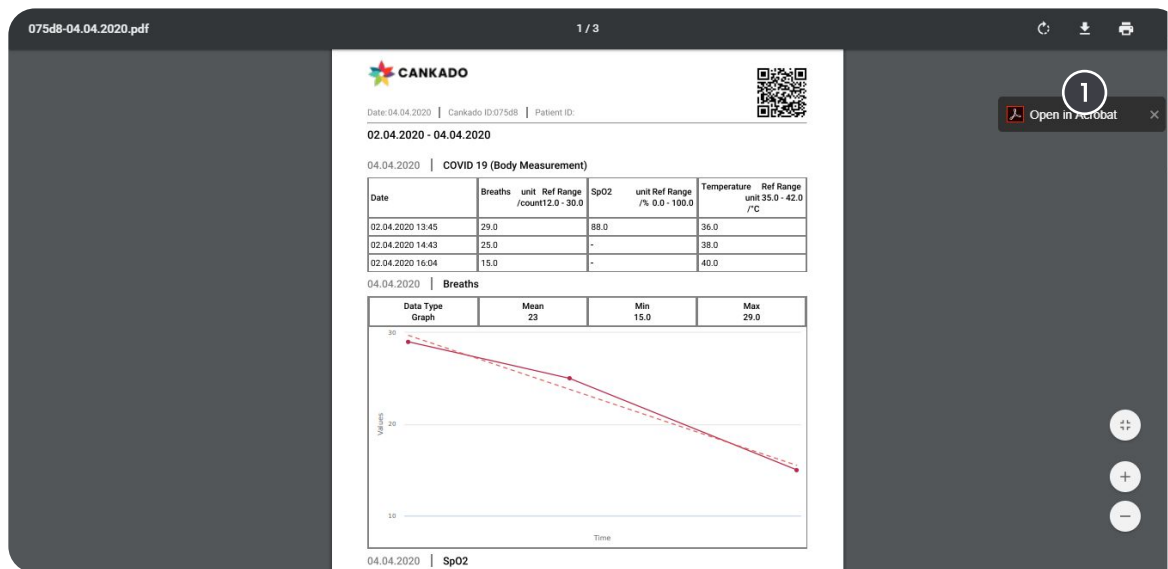

After selecting [6] in Cockpit Overview as depicted on page 3 the corresponding patient report will open as pdf file.

1) Click here to download the file.

7: Export patient list

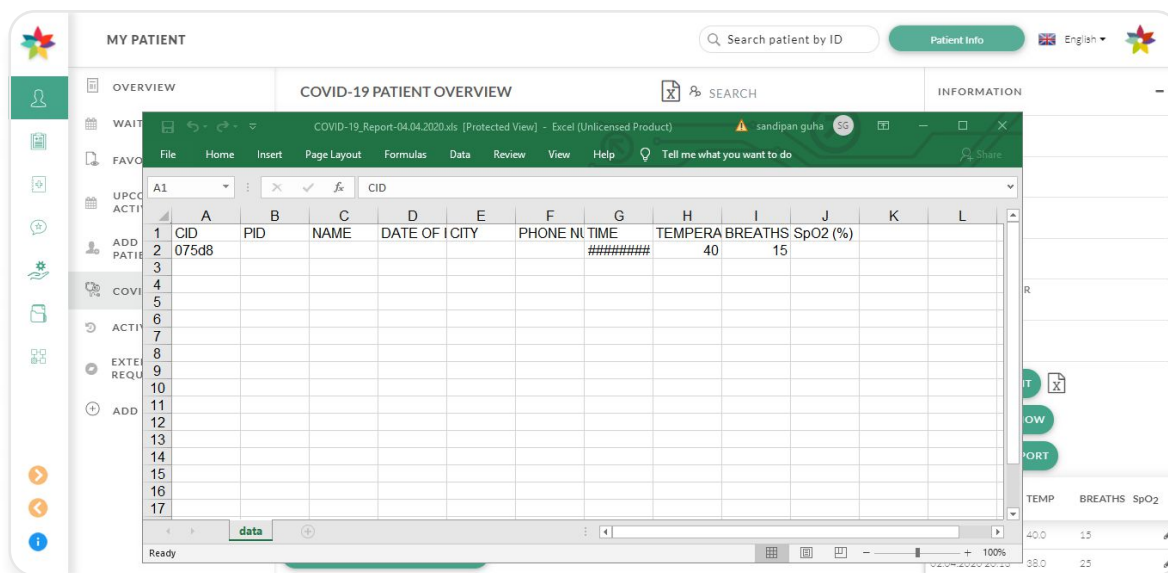

After selecting [7] in Cockpit Overview as depicted on page 3 your patient list will be exported as xls file.


8: Go to patient or view patient information

The screenshot displays the 'MY PATIENT' interface. On the left is a sidebar with navigation options: OVERVIEW, WAITING ROOM, FAVOURITES, 075D8, UPCOMING ACTIVITIES, ADD EXISTING PATIENT, COVID-19, ACTIVITY HISTORY, EXTENSION REQUESTS, and ADD PATIENT. The main area is titled 'PATIENT INFORMATION (075d8)' and contains a 'Personal Information' section. This section includes fields for PATIENT ID (PID), DATE OF BIRTH, GENDER, FIRST NAME, LAST NAME, EMAIL, PHONE NUMBER, STREET, CITY, POSTAL CODE, and COUNTRY. There are two buttons: 'ASSIST' (labeled 4) and 'WELCOME SHEET' (labeled 2). A '1' is placed over the 'Select date' field. A '3' is placed over the 'COVID-19' option in the sidebar. The right sidebar shows 'BASIC SETTINGS' with checkboxes for 'Adherence, only if treatment is assigned', 'Neuropathy Pain Scale (0-10)', and 'Global Health Status (0-100)'. Below this is 'PATIENT SETTINGS OVERVIEW' with a checkbox for 'Enable patient settings' and a 'SAVE' button.

After selecting [8] in Cockpit Overview as depicted on page 3 the corresponding Patient Information page will open.

- 1) Here you can add or edit patient ID, date of birth, name and address.
- 2) Here you can print out a "Welcome Sheet" document for your patient containing general information and patient's login data.
- 3) Here you can go back to "Covid-19 Cockpit Overview".
- 4) Here you can open up "Assist Mode". Activating this mode will switch you to this patient's account, as if the patient himself had logged in. To end this mode select the little Cankado symbol in the upper right corner and click "Logout".


How to connect to a patient

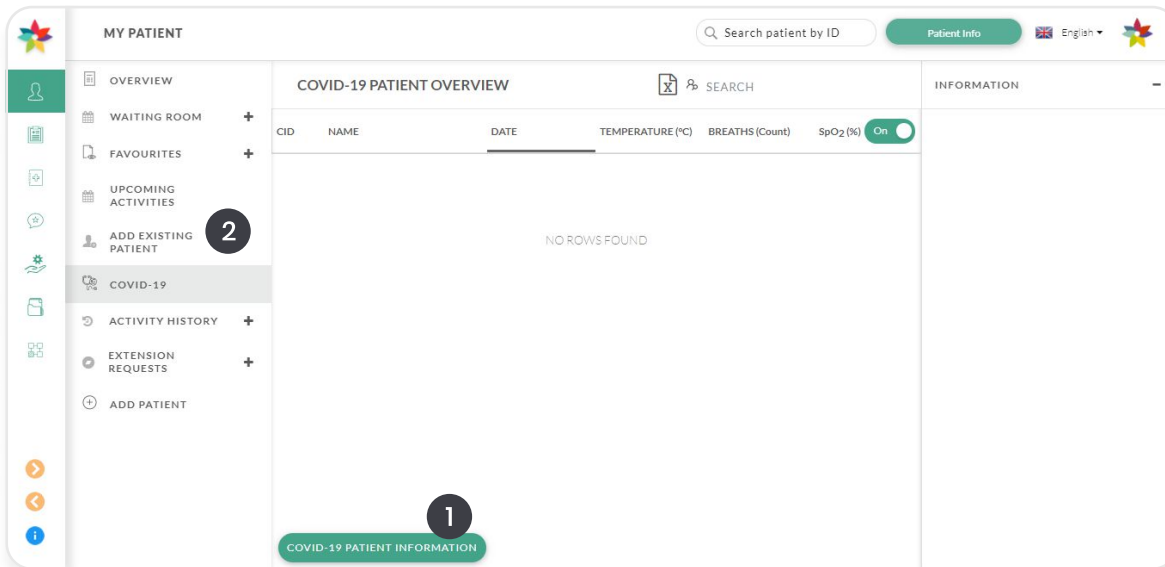

There are several ways to connect to a patient.

- 1) Use this option if your patient is a new user or an existing user who has not activated Covid-19 Modul.
- 2) Use this option if your patient is an existing user and already activated Covid-19 Modul.

1: Connect patient via QR code provided

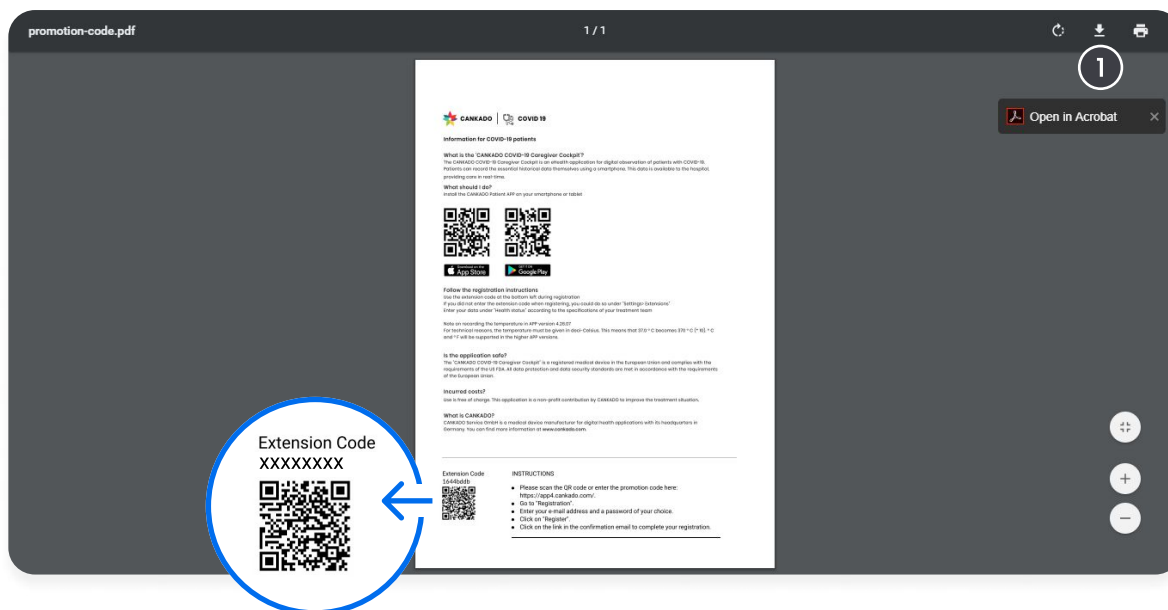

After selecting [1] in Cockpit Overview as depicted on page 7 the “COVID-19 Patient Information” document opens. Scanning the bottom left QR-Code will enroll your patient in the Covid-19 Modul and add the patient to your department in one step.

- 1) Click here to download the document (to print it out or send it to your patient).

2: Add existing patient


The screenshot shows the 'MY PATIENT' interface. On the left is a sidebar with icons for Overview, Waiting Room, Favourites, Upcoming Activities, Add Existing Patient (selected), COVID-19, Activity History, Extension Requests, and Add Patient. The main area is titled 'ADD EXISTING PATIENT' and contains a search bar 'Search patient by ID', a 'Patient Info' button, and a language dropdown set to 'English'. Below the search bar is a text input field labeled 'CODE TO ADD PATIENT' with a red asterisk, a green 'SAVE' button, and a circled '1' indicating the input field. On the right is an 'INFO' panel with a yellow instruction box: 'Enter the activation code from the physician invitation sheet to enroll the patient to the current department.'

After selecting [2] in Cockpit Overview as depicted on page 7 this windows opens.

- 1) Type in the Invitation Code you received from your patient (for information on where your patient can find this code see page 13 of our "Patient Instruction Manual").
- 2) Click "Save" to continue.

How to register


Visit <https://cankado.com/>

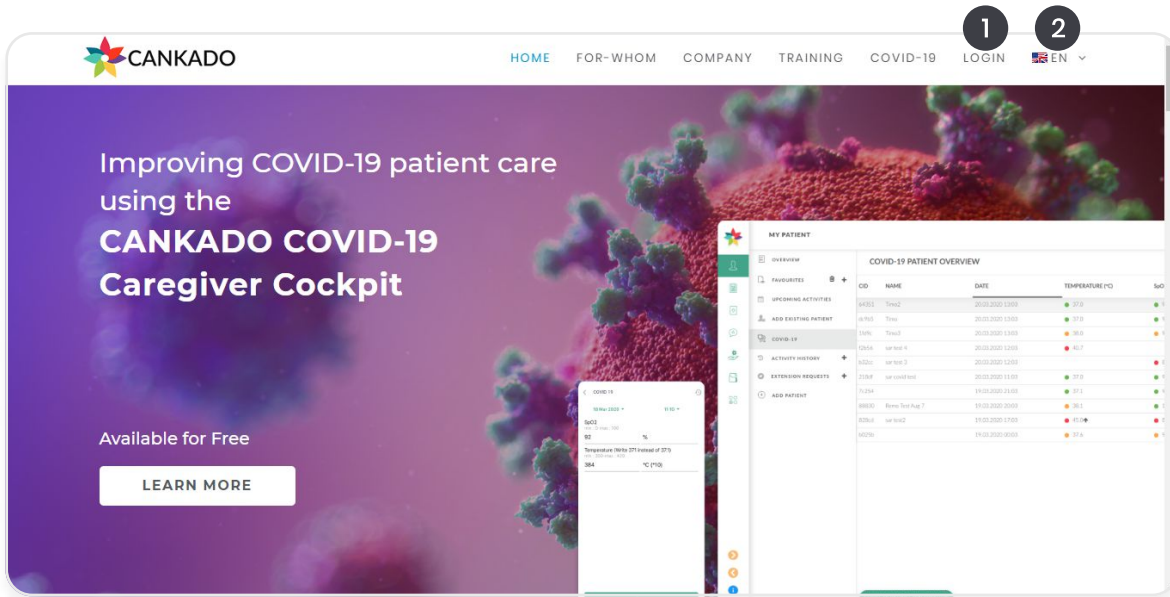

Open cankado.com in your web browser.

- 1) Click here to login.
- 2) Click here to switch language.

Select **“Register”** from top

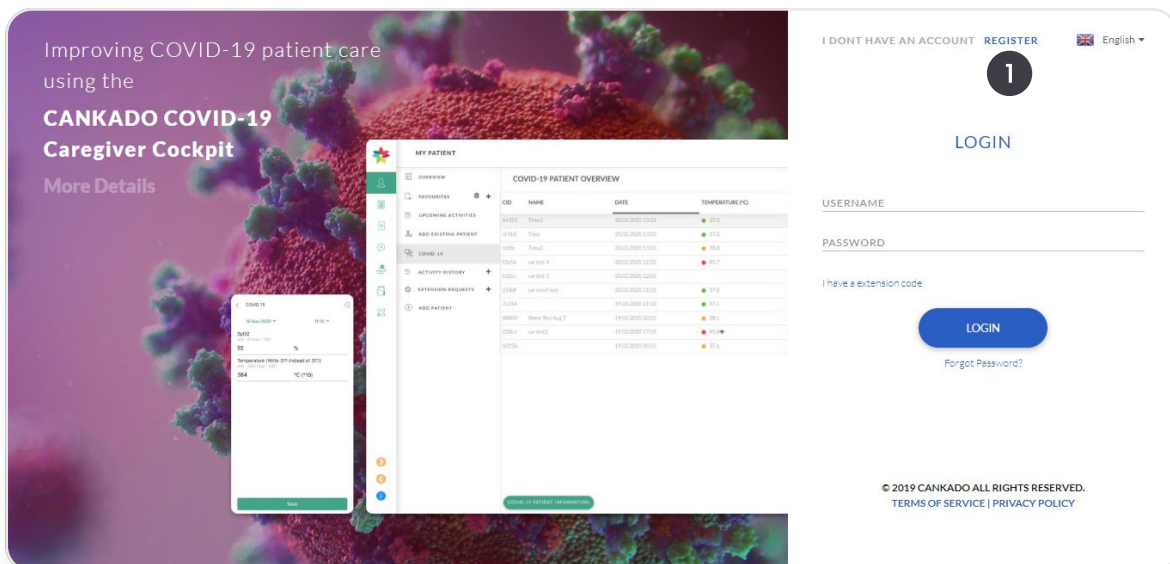

- 1) Select **“Register”** on the top of the page.

Select doctor registration

The screenshot shows the CANKADO website's doctor registration process. On the left, a blue banner reads "THANK YOU FOR YOUR INTEREST IN CANKADO" and "REGISTRATION FOR PHYSICIAN". On the right, the "DOCTOR REGISTRATION" form is displayed with four numbered steps: 1) "ALREADY A MEMBER LOGIN" at the top right; 2) "EMAIL" input field; 3) "PASSWORD" and "RETYPE PASSWORD" input fields; 4) "VERIFICATION" dropdown menu with "Efn verification" selected. Below the dropdown is a "REGISTER" button. At the bottom, there is a disclaimer: "BY CLICKING REGISTER YOU AGREE TO THE TERMS AND CONDITIONS AND OF THE PRIVACY POLICY." and a copyright notice: "© 2019 CANKADO ALL RIGHTS RESERVED."

- 1) Select "Doctor" as user type.
- 2) Type in a valid email address. Choose a strong password and re-enter it to verify.
- 3) Physicians using CANKADO must be verified. Because of the urgency of the situation, we have simplified our validation process. Register with an email address listed on your hospital homepage and select "Manual Verification". We will activate your account as soon as possible.
- 4) Select "Register" to continue.

Welcome to homescreen & add a canter

The screenshot shows the CANKADO "WELCOME TO CANKADO" screen. It features the CANKADO logo at the top left and a language selector at the top right. The main heading is "WELCOME TO CANKADO" with the subtext "What do you want to do?". Below this, it says "SEEMS LIKE YOU HAVE NO DEPARTMENTS". A prominent orange button labeled "GO TO CENTER MANAGEMENT" is centered, with a numbered step "1" next to it. At the bottom left, there is a small copyright notice: "© 2020 All rights reserved."

First you have to add a center and a department to your new physician account.

- 1) Click on "Go to Center Managment".

How to add a center

The screenshot shows the 'CENTERS' section of the Cankado interface. On the left, there is a sidebar with 'OVERVIEW' and 'ADD NEW CENTER' options. The 'ADD NEW CENTER' form is displayed with two input fields: 'CENTER NAME' (marked with a red asterisk and a circled '1') and 'ADDRESS'. A green 'ADD' button is located at the bottom right of the form, marked with a circled '2'. The top of the interface includes a search bar for patients by ID, a language dropdown set to 'English', and a user profile icon.

- 1) Enter your center's name here.
- 2) Continue with "Add".


How to create a department once a center is created

The screenshot shows the 'CENTERS' overview table and the 'CENTER INFORMATION' panel. The table has columns for 'CENTER NAME', 'TYPE', 'DEPARTMENTS', 'OWNER', and 'ACTION'. A row is visible for 'COVID_Track_010' with 'General' type, '0' departments, and 'Myself' as the owner. A 'DELETE' button is in the 'ACTION' column. The 'CENTER INFORMATION' panel on the right shows details for the selected center, including 'CENTER NAME', 'CENTER TYPE', 'OWNER', 'PHONE NUMBER', and 'HOME PAGE'. A green 'MORE INFO' button is at the bottom of the panel, marked with a circled '1'. The top of the interface includes a search bar for patients by ID, a language dropdown set to 'English', and a user profile icon.

CENTER NAME	TYPE	DEPARTMENTS	OWNER	ACTION
COVID_Track_010	General	0	Myself	<button>DELETE</button>

- 1) To add a department to your new center click on "More Info".


Create a department

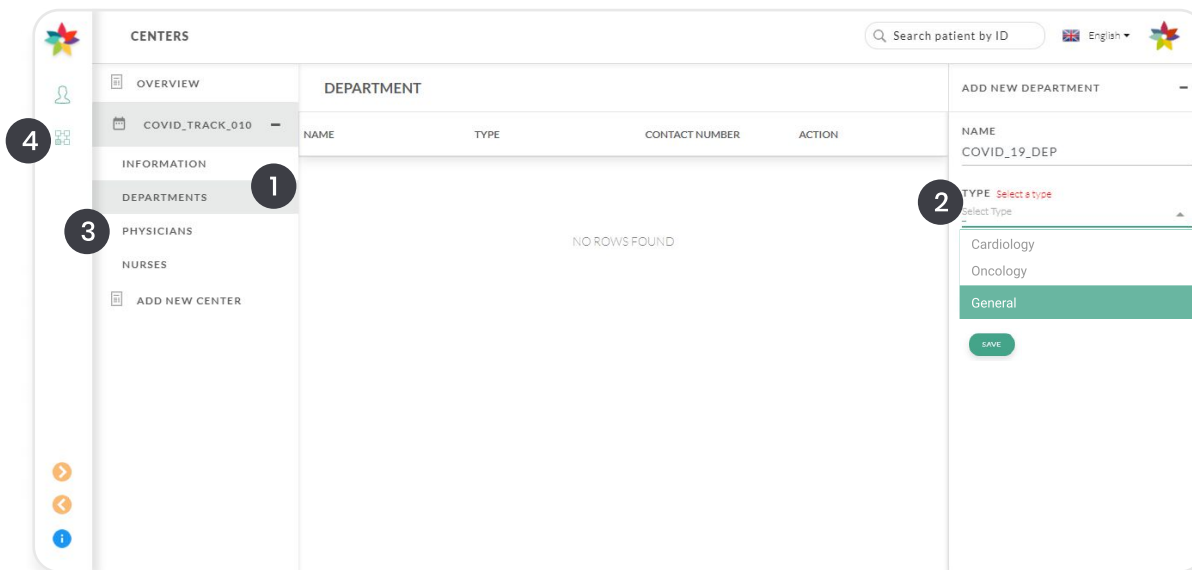

- 1) Click on "Departments".
- 2) Enter a name and choose "General" as type. Click on "Save".
- 3) If you want to invite other physicians to your department click on "Physicians" and continue with next step.
- 4) If you don't want to invite other physicians at this point (you can still do this later on) click on "Departments".

"Important Note: The COVID-19 Module ist only available in department types general and oncology. Do not use any of the others."


Invite your fellow phisicians

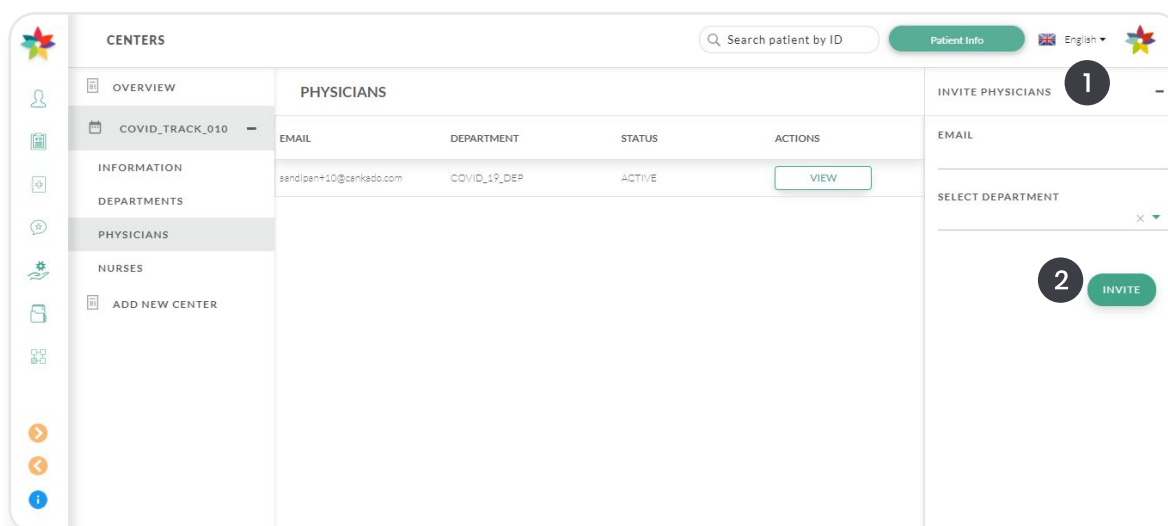

Continue here if you want to invite other physicians to your department and have selected [3] in the top picture of this page.

- 1) Enter email address of the physician you want to invite.
- 2) Select "Invite".
- 3) Click on "Departments" to continue.

Go back to home screen and explore COVID-19 cockpit

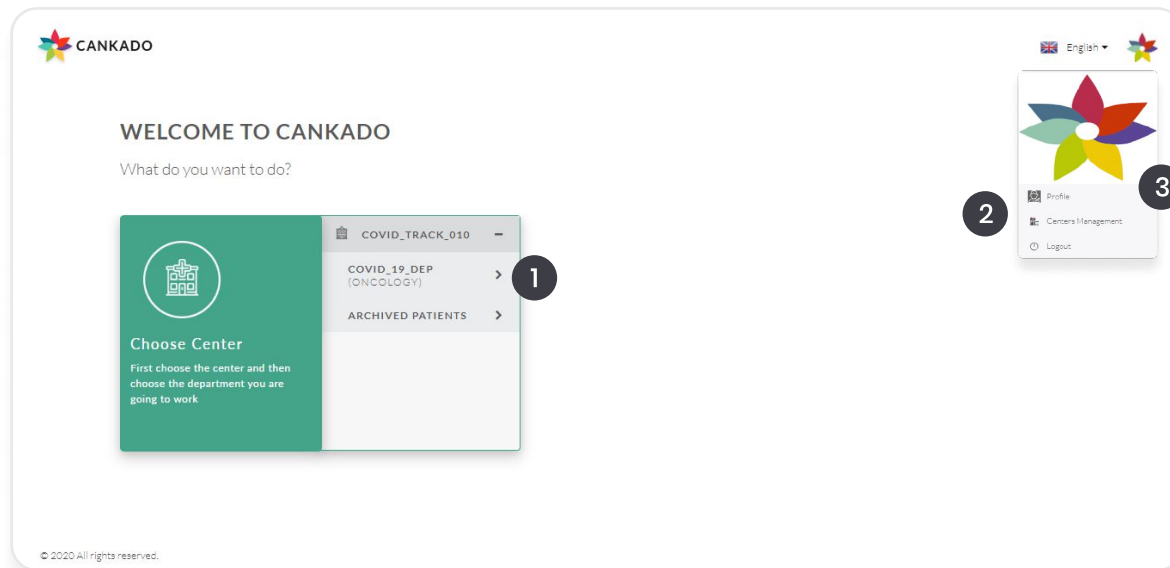

Selecting “Departments” from the left side menu brings you to this screen. Here you can see an overview of your centers and departments.

- 1) Click on “Chosse Center” and then select your newly created center and department to enter.
- 2) To go back to “Center Management” click here.
- 3) Click on “Profile” to set your new center and department as your default department. This step is needed for adding a patient to your department as depicted on page 19.

Go to Profile and set your department as Default Department

To add a patient to your department using the patient’s QR code as depicted on page 19 you need to set a default department.

- 1) After selecting “Profile” as shown in the previous picture, click on “Departments”.
- 2) Select your center and department here.
- 3) Click “Save” to confirm.



Using the HCP Pro App on a smartphone

App overview

HCP Pro home

Here you can see the home screen of HCP Pro App.

- 1) Click on "COVID-19 Monitor" to see an overview of all your COVID-19 patients.
- 2) Click here to see an overview of all your patients in the selected department.

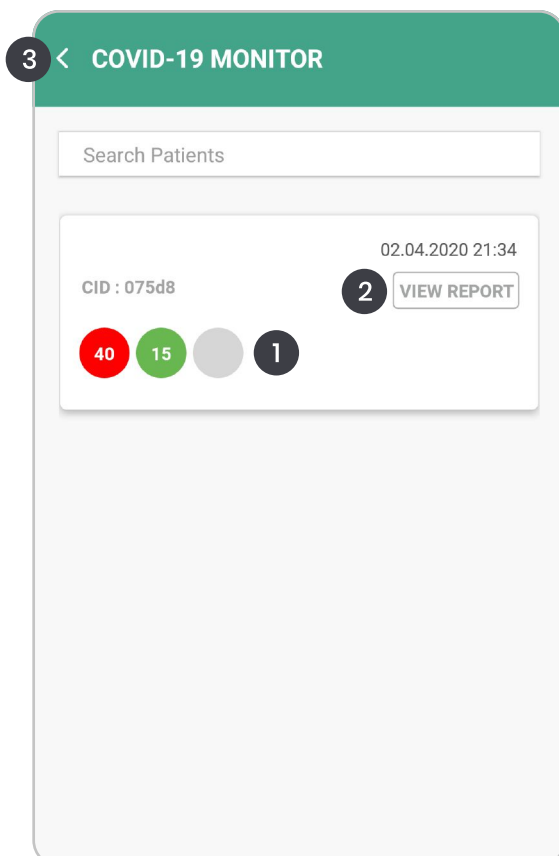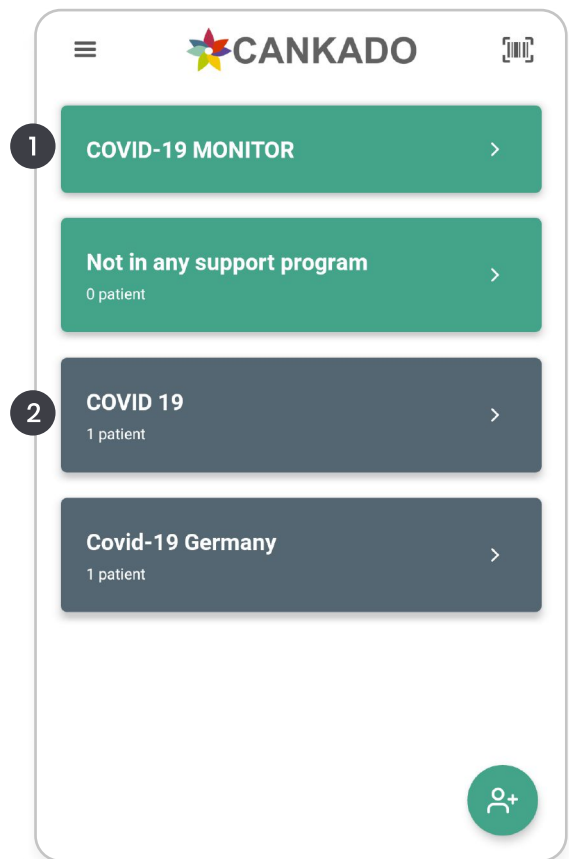

COVID-19 monitor

Colour code	Temperature	Breaths	SpO2
green	< 37.5 °C	< 19	> 94 %
Orange	37.5 - 38.5 °C	19 - 26	90 - 94 %
red	> 38.5 °C	> 26	< 90 %

1) Here you can see the values last documented for the corresponding patient. For information regarding color code see table depicted above.

Click on it to see all documented values for this patient.


2) To see patient's report click on "View Report".

3) Click here to go back to the previous screen.

Generate PDF report


Selecting [2] in the last step will open up the Patient Report document.

- 1) Here you can see the documented values as report graph.
- 2) Select "Share" to share this document (e.g. per email).
- 3) Click here to go back to the previous menu.

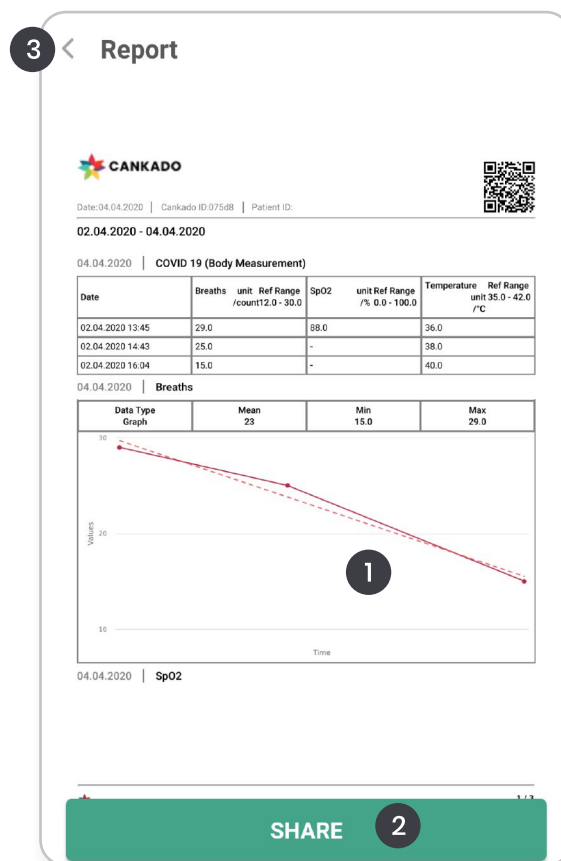

Share report

After selecting [2] in the previous step you can choose in which way you want to share the report.

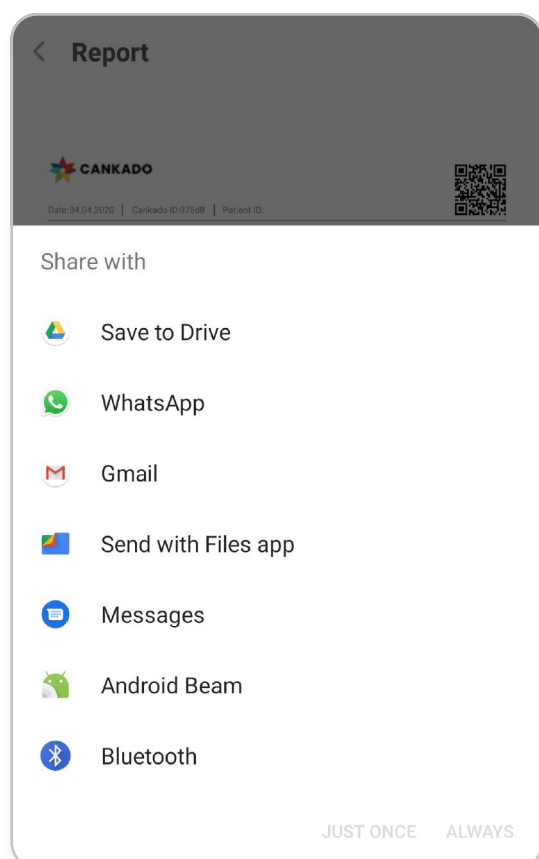

View Patient data


Selecting [1] in Covid-19 Monitor as depicted on page 14 opens up a table containing all documented data.

- 1) Here you can see date of documentation.
- 2) Here you can open up graph view.
- 3) Click here to go back to the previous screen.

3 < 075d8

2

Time	Temp.	Breaths	SpO2
04.04.2020			
17:28	40	15	-
03.04.2020			
15:20	41	18	-
02.04.2020 1			
21:34	40	15	-
20:13	38	25	-
19:15	36	29	88

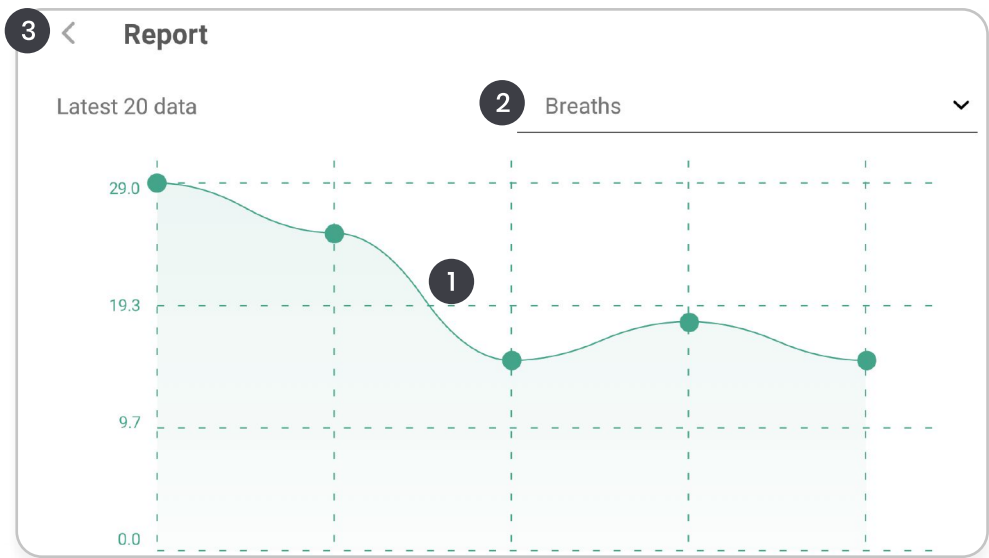



View graph

- 1) Here you can see the report graph for the selected value.
- 2) Click here to select a different value.
- 3) Click here to go back to the previous screen.

List of all patients

Selecting your department [2] on the home screen as depicted on page 14 opens up this view. Here you can see an overview of all the patients added to your department whether these patients have an active COVID-19 module or not.

- 1) Click here to select corresponding patient.
- 2) Type CID or patient name here to search for patient.
- 3) Click here to go back to previous screen.

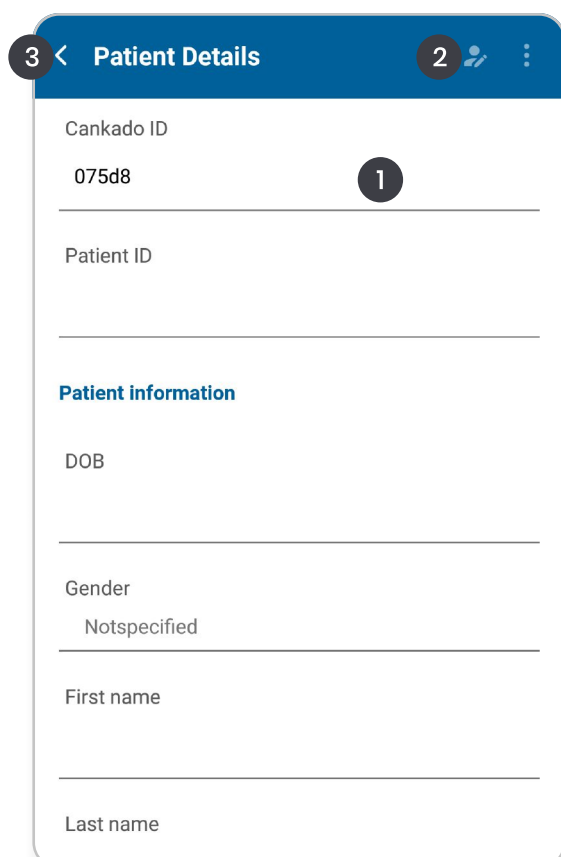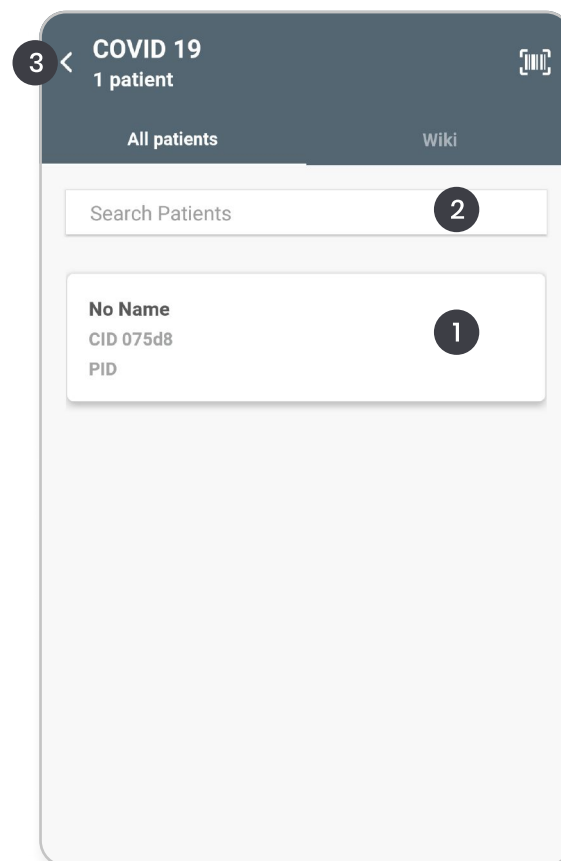


View patient details

After selecting [1] Patient Details will open.


- 1) Here you can see patient's CID.
- 2) Click here to edit Patient Details.
- 3) Click here to go back to previous screen.

Push notification

Whenever a patient documents values within orange or red zone for the first time you will get a push notification informing you about it.

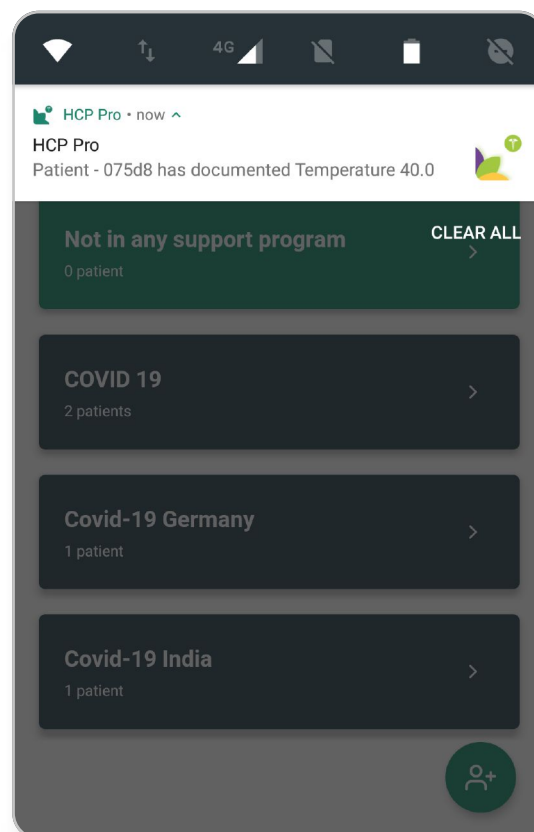


How to connect to a patient

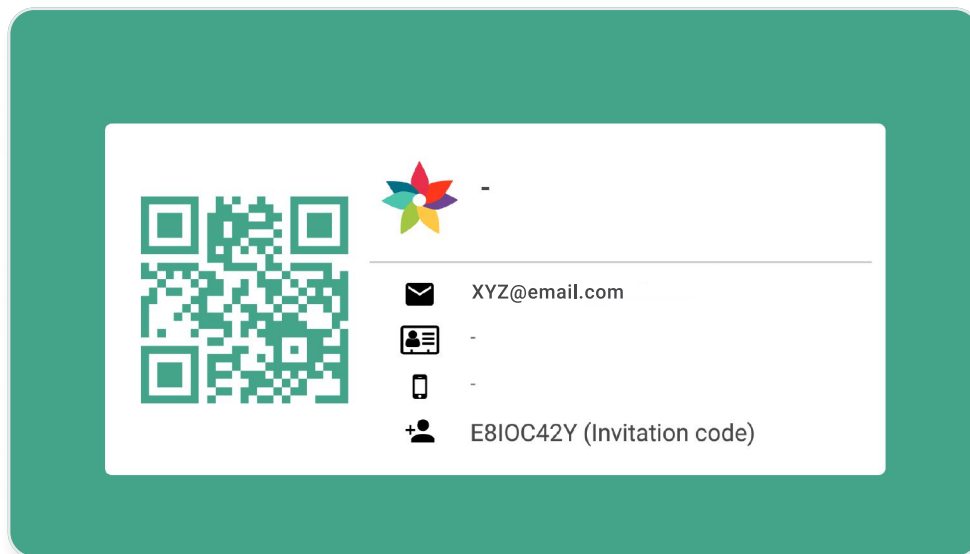

Patient QR code

To connect to a patient you can scan the patients QR code with your own HCP Pro App. For information on how your patient can find this code see page 13 of our "Patient Instruction Manual".

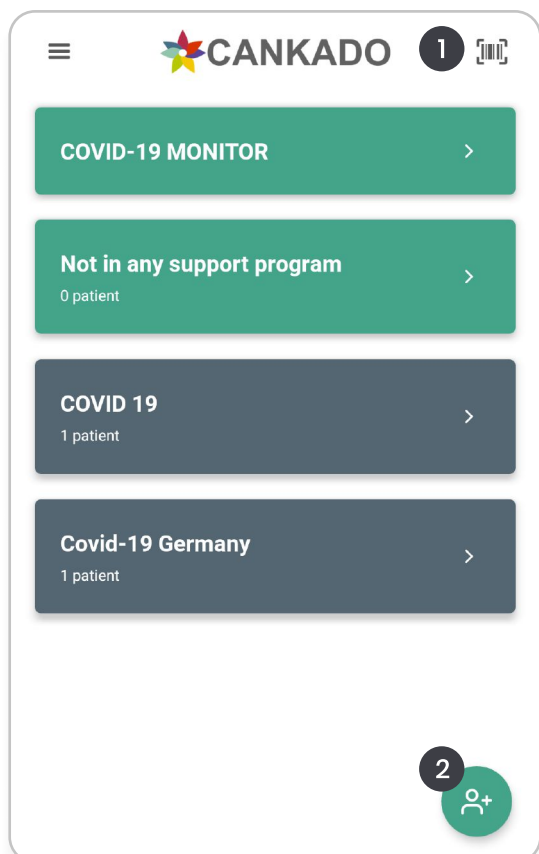


Select QR code scanner from home

1. Click here to scan a QR code and connect with an existing patient.
2. Alternatively you can click here to add a new patient.

Scan QR code

If you selected [1] in the previous step you can now scan the patient's QR code.

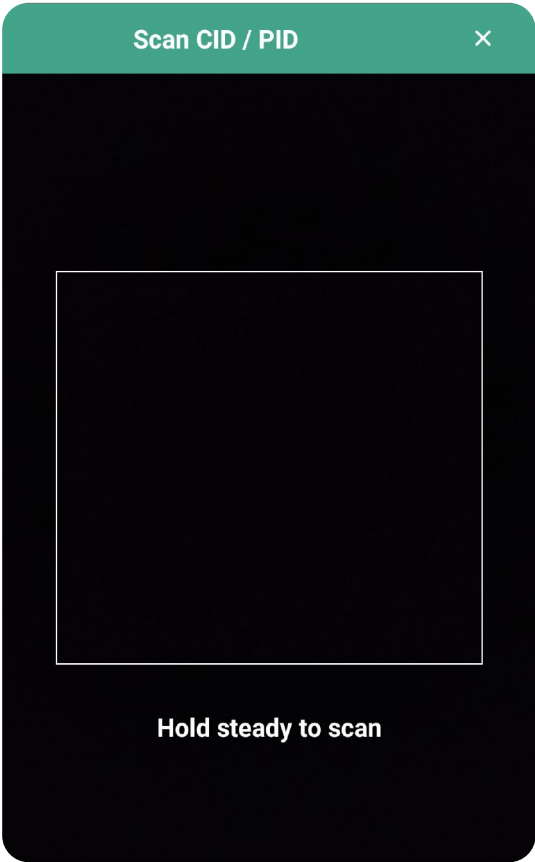

×

Add new patient

Patient ID

Enter PID

1

2

Date Of Birth

Select Date of Birth

Gender

Choose an option.

▼

First name

Enter first name

Last name

Enter last name

Email ID

Add as new patient

Go to add new patient

If you selected [2] in the previous step you will see this view.

- 1) You can either add a completey new patient,
- 2) Or you can scan the QR code of an existing patient.

Scan QR code

You can now scan the patient's QR code.

Add new patient

Patient ID

Enter PID

Scan QR Code / Barcode

1


x

Email ID

Add as new patient

How to register

Visit <https://cankado.com/>

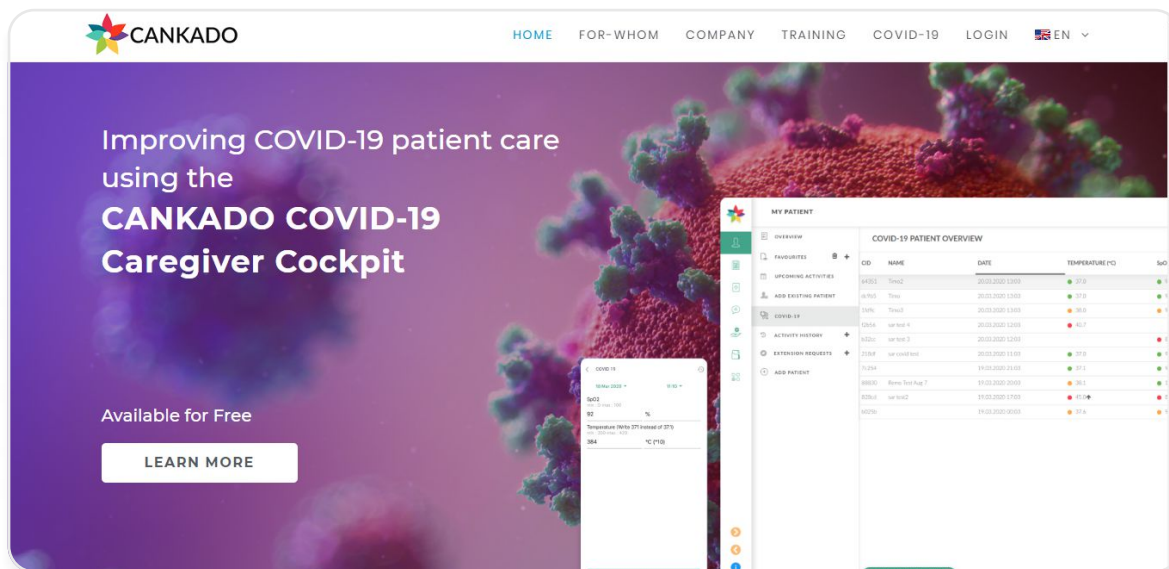

Go to cankado.com and follow the same steps shown in page [8].



CANKADO

Health
made in Germany

Contact Us

+49 (0) 221 429 153 00
covid19@cankado.com

OUR OFFICES

Munich

Am Strassland 6
85551 Kirchheim b.M
Germany

Buenos Aires

Bahía Blanca 2390
Buenos Aires
Argentina

Trivandrum

Technopark
Campus, Amstor
Kerala 695582
India

Cologne

Eupener Str.165
50933 Cologne
Germany

Indianapolis

Crosspoint Plaza One
10475 Crosspoint Blvd
Indianapolis, IN 46256
USA